# CHANGE X ARM

USER MANUAL

- › English
- › Norsk
- › Dansk
- › Svensk
- › Suomi


EXPAIN®

## Contents

- 1 Components
- 2 Introduction
- 3 Before use
- 4 How to mount the device
- 5 How to unmount the device
- 6 Usage
- 7 Usage with app
- 8 Biofeedback and habits
- 9 Battery and charging
- 10 Maintenance and storage
- 11 Technical specifications
- 12 Safety information

## Components

1. Button
2. Light
3. USB connector
4. Adjustable arm strap
5. USB cable
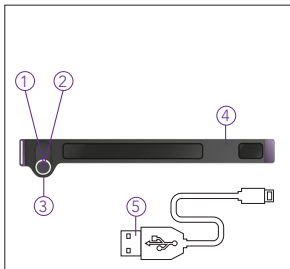
## Introduction

Congratulations on your new **EXPAIN Change Arm!** This is a unique device that allows you to monitor the muscle activity in the forearm muscles in a new way. **EXPAIN Change Arm** provides Biofeedback to help you change habits that may cause pain. Biofeedback is a learning process where observation of the signals from your own body, leads to a changed attitude. The device vibrates when you have strained muscle activity pattern. This allows you to change your bad habits, which again may eliminate the cause of your pain.

**EXPAIN Change Arm** can be connected with a smartphone via Bluetooth Smart, or it can be used without connection to phone. If you wish to keep track of your historical muscle activity pattern, you must

download **EXPAIN app** (iOS or Android). See separate section or [www.expain.com](http://www.expain.com) for more information.

The unit is wireless, powered with a rechargeable battery and is placed on the painful muscle in the forearm, and where the forearm is at its thickest.

### **Before Use**


**EXPAIN Change Arm** is made for personal use. Read instructions carefully. The unit must be charged before the first use. See section for battery and charging for more information on how the device is charged. *(NB: The device should not be used at night, and the unit should not be on the body continuously for over 16 hours. If there are signs of allergic reactions or rash, immediately stop using the device. Use the device with caution).*

## How to mount the device

1. Loosen the lock function on the armband and pull it over the hand and up onto the forearm.
2. Make sure the rubber fields inside of the belt are over the most painful muscle and across the muscle fibers. Attach the product on the side you have the most pain, for most people this is the dominant side (right for right-handers and left for left-handers).  
If you wish to measure both arms simultaneously we suggest one

3. extra device.  
Tighten the belt so hard that it is snug but does not restrict movement or blood flow in the forearm.

Let the device be on your body for the time period you want biofeedback.
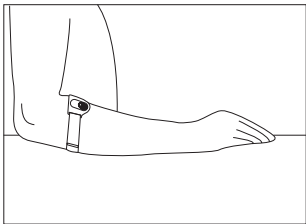
## How to unmount your device

1. Take the unit of the body by loosening the belt tensioning function and take it off like a bracelet.
2. Recharge the unit to ensure that it is ready for the next use.

## Usage

### **EXPAIN Change Arm** is

turned on by pressing the button on top of the unit.

When using the device, it will vibrate when it recognizes strained muscle activity patterns.

Read more in the section of Biofeedback and habits.

### Usage with app

You can get more from the **EXPAIN Change Arm** device if you connect it to a smartphone.

The unit will then automatically transfer data to your smartphone when it is within reach and you can follow the historical development in your app. You can, for example do ergonomic changes by looking at your muscle activity.

1. **EXPAIN Change Arm** is connectable whenever the device is turned on.

It is not possible to connect the device to multiple phones simultaneously. If you want to pair the device with another phone, you must remove the device in the settings in your app / phone, and restore factory settings if necessary, see the Battery and

charging section.

The **EXPAIN app** can be downloaded from Apple App Store or Google Play Store. Follow the instructions in the app. The app has a built in pain diary that allows you to compare your muscle activity patterns with your pain. In this way you can see when you've had the most strained muscle activity and if this has a relationship with your pain. It is recommended to register your pain regularly, and to use the device at least for a few days to get the most out of **EXPAIN Change Arm**.

## Light and vibration

The table below shows the meaning of various lights and vibrations.

Function	Light/Vibration
Connected Bluetooth Smart	Blue flash
Low battery	Red flash
On, without Bluetooth connection	Purple flash
Low battery	Red flash
Fully charged	Purple
Notice sustained muscle activity	Short vibrations
Low battery	Long vibrations

**Biofeedback and habits**

Habit change happens when the device is used often and you take breaks or do some exercises every time the device vibrates. Research shows that biofeedback, which is a learning process, helps to improve health and overall performance by observing the signals generated by one's own body.

The device will vibrate when it recognizes strained muscle activity patterns (Sustained muscle activity). When the device vibrates, try to take a short break from work, do some arm exercises or lower your hands so that the muscles in the forearm region relax. A good relaxation exercise is to breathe calmly in with the nose and then releasing the

breath out of your mouth while you let your arms and shoulders hang down by the side of the body. Repeat this exercise 2-5 times.

### **Battery and charging**

The device has a battery life of approximately 8 hours of normal use. The sensor flashes red to indicate low battery. In addition, the device vibrates Long vibrations when the battery is critically low.

**EXPAIN Change Arm** is charged with supplied or corresponding micro USB cable. Unmount the device from your arm and charge it completely before next use. The lights is red during charging, and purple when it is fully charged.

Restoring factory settings. If you can't connect your phone the device can be reset to its original

state. This is done as follows: Turn the device on and set it to charge. Press the button until the device starts flashing red. When the device lights red and vibrates, the factory settings are restored.

### **Maintenance and Storage**

Clean the device with a dry or slightly damp cloth.

Make sure the device does not get wet. It is not recommended to use the device in humid environments or if you plan to do activities that will cause sweating.

The Arm-strap can be washed in water, but only after that the electronic unit has been removed.

The unit should be transported, stored and operated in temperatures between 10°C and 40°C, at a relative humidity

between 30% and 85% and with an atmospheric pressure between 500 and 1060 hPa. If the device does not work correctly, contact support at [www.expain.com](http://www.expain.com)

Tehnical specifications	
Technology	Surface Electromyography (sEMG)
Power	LiPo rechargeable battery
Bluetooth	Bluetooth 4.2 Smart
Size (DxBxH)	65 mm x 360 mm x 15 mm
Weight	45 g
Safety	IEC 62133 EN 60950
EMC	EN 301 489-17 EN 300 328 EN 301 489-1
Safety	RoHS, REACH Primary Skin Irritation Test
Made in	Norway

## **Safety information**

### Applications

For Biofeedback against strained (sustained) muscle activity in the forearm muscles.

### **Safety precautions**

Do not use this product if you have very sensitive skin.

### ***Warning!***

*Do not use the device on other locations on the body than those specified in this manual.*

**Expain AS**  
Karl Johans gate 13,  
0154 Oslo, Norway  
[www.expain.com](http://www.expain.com)

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