EXPAIN® RELAX MINI

USER MANUAL

- › English
- › Norsk
- › Dansk
- › Svensk
- › Suomi

EXPAIN®

Contents

- 1 Introduction
- 2 Components
- 3 Before use
- 4 Wheel massage
- 5 Vibration massage
- 6 TENS function
- 7 Battery and charging
- 8 Maintenance and storage
- 9 Technical specifications
- 10 Safety information

Introduction

Congratulations on your new **EXPAIN Relax Mini**. This is a unique and innovative 3 in 1 unit, with spot massage, vibration massage and TENS treatment.

This innovative device is safe, without medications, easy to use, discreet in the purse, comfortable to use - and it gives you the opportunity to live an active life.

EXPAIN Relax Mini consists of 4 innovative massage heads that have different functions.

Components

PARTS


- 1 Wheel Massage
- 2 Vibration Massage
- 3 TENS points
- 4 LED and micro-USB connector

CHARGING

- 5 USB Charge Cable
- 6 Control ring

CONTROL

- 7 Control Ring rotated to set mini massage to TENS (Low, medium and high intensity) vibration or off.
- 8 Interface with LED indicates charging status, TENS and vibration.

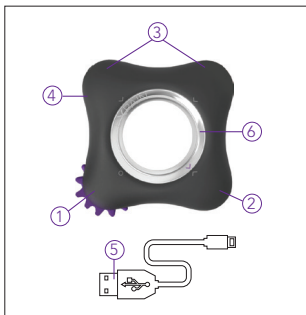Before Use

EXPAIN Relax Mini is designed to be used as needed. It is intended for a massage or a quick electric stimulation.

Read the manual and safety precautions before use.

EXPAIN Relax Mini can be used on battery or by power cord.

Charge your device fully before first time use. Read more in section “Battery and Charging.”

Wheel Massage

Manual spot wheel massage that you can roll back and forth on your muscles.

1. The arrow in the control ring should be pointing to the O symbol, which is towards the wheel.
2. Place expain Relax Mini on your muscles.
3. Scroll back and forth in the direction of the muscle fibers.

Vibration Massage

The Vibration Massage head gives shock massage by vibration of 6500 RPM (rotations per minute). This provides optimum power per area and creates flow to the muscles.

1. Turn on vibration massage by turning the arrow control ring to the massage head having vibration therein (2).
2. LED lights purple when vibration is in progress.
3. Turn the control ring to O symbol to stop the vibration massage.

TENS function

TENS does trans electrical nerve stimulation. TENS technology has been used by therapists in several decades to relieve pain.

1. Turn the control ring so that the arrow on the ring points to the lowest intensity of TENS.
LED lights blue when TENS is on.
2. Place the two massage heads ,which have black ends, on the skin in the area where you have

aching muscles.

3. Now you can turn the control ring to increase the intensity of TENS activation. Don't use higher intensity if painful.

The effect can be increased by wetting the black electrode areas with water.

Warning!

TENS should not be used if you have a pacemaker, are pregnant or have heart problems.

Ask your doctor if you are unsure.

Do not apply TENS in the face, temples or crotch.

Turn off the TENS function by turning the control ring so that the arrow points to the O symbol. TENS turns automatically off after 15 minutes.

It is recommended to have a minimum of a 30 minute break between two 15

minute TENS sessions. EXPAIN Relax Mini includes a USB charging cable and Power outlet. The device has a battery life of approximately 1 hour with normal use.

The device shows low battery with a red flashing light. Charge the unit using the supplied cable or continue to use the device with charging cable connected. When charging the light is red. When the device is fully charged the Light is purple, and it's ready for battery use again.

Maintenance and Storage

Clean the device with a dry or slightly damp cloth.

Keep the device out of reach of children and animals. **EXPAIN Relax Mini** is not waterproof. Do not immerse **EXPAIN Relax Mini** in liquids. The unit should be transported, stored and operated

within temperatures of 0 ° C and 40 ° C, at a relative humidity between 30% and 85% and with an atmospheric pressure at 500-1060 hPa. If the device is not working properly, stop all use and seek support on www.expain.com

Technical specifications

Massage Technology	Wheels and vibration massage
TENS-Technology	Automatic modulating frequency program
Power	LiPo battery or cable
Size (DxBxH)	12 cm x 12 cm x 4 cm
Weight	0,3 kg
Safety	EN 60335 UN38.3 RoHS REACH
Made in	China

Security info

Areas of use

For use on muscles.

Safety rules

Do not use this product if you have open wounds or very sensitive skin that can't tolerate massage or TENS. Do not use this product if you have muscles that can't withstand normal manual massage.

Warning!

Do not use the device in the face, temples or crotch.

Do not use TENS function if you have a pacemaker, implanted defibrillator or heart problems. Consult your doctor if you are unsure.

Do not use expain Relax Mini on stomach if you are pregnant.

TENS activation for other activities may increase the risk of damage. Do not use the device when bathing, showering, sleeping, running, controlling machines or perform other activities that could be disturbed by TENS activation.

*Using TENS function Near instruments (as eg ECG equipment, **EXPAIN Change Neck** or **EXPAIN Change Arm**, ECG alarms) may cause the Equipment to stop working.*

*Do not use TENS function near electrical measuring instruments. The efficacy and safety using TENS function in children, pregnant women, or the head have not been investigated. Do not use TENS function on children or allow children to use **EXPAIN Relax Mini**.*

Do not use **EXPAIN Relax Mini** if you are pregnant, or suspect that you are pregnant, unless your therapist recommends it. Do not use **EXPAIN Relax Mini** on the head.

If you have or suspect you have epilepsy or heart disease, TENS activation may have unexpected effects. Always ask a doctor before use of TENS activation if you have epilepsy or heart failure.

Expain AS
Karl Johans gate 13
0154 Oslo, Norway
www.expain.com

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