EXPAIN® RELAX NECK AND SHOULDERS

MANUAL

- › English
- › Norsk
- › Dansk
- › Svensk
- › Suomi

EXPAIN®

Contents

- 1 Components
- 2 Introduction
- 3 Before use
- 4 Massage function
- 5 Heating function
- 6 Batteries and charging
- 7 Maintenance and storage
- 8 Technical specifications
- 9 Safety information


Components

PARTS

1. Neck massager
2. Button On / Off / rotation
3. User interface
 - a. LED
 - b. Micro-USB connector for charging
4. USB charge cable
5. Handle

CONTROL

Button	Function	Light
Hold for 2 seconds	On/Off	Purple
Short press	Soft massage	Purple
Short press	Hard massage	Blue
	Low battery	Red flash
	Charging	Red
	Charged	Purple

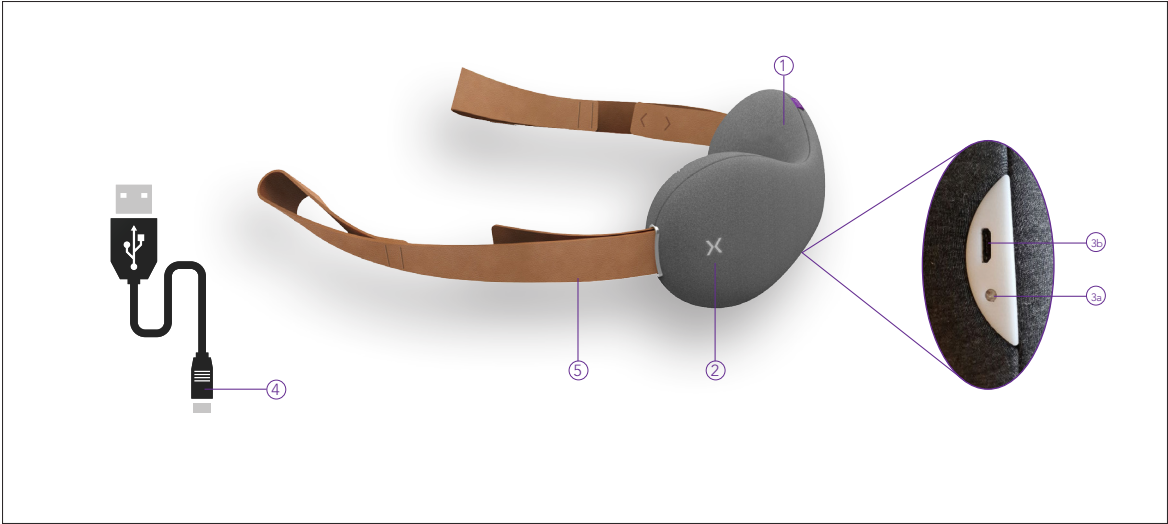Introduction

Congratulations on your new **EXPAIN Relax Neck & Shoulders**. This is a unique and battery-powered device for the neck and shoulders. **EXPAIN Relax Neck & Shoulders** has hard and soft massage and active heat. This innovative device is safe to use, medication free, easy to use, discreet in the couch, comfortable - and above all, it gives you the opportunity to live an active life.

Before Use

EXPAIN Relax Neck & Shoulders is designed for personal use. Read the instructions and safety precautions

before use.

You can use **EXPAIN Relax Neck & Shoulders** on battery or with a cable to the wall. Charge the device for 2 hours before the first time you start using it. Place the transformer plug in the socket. Insert the USB cable in the transformer plug, and then insert it in the **EXPAIN Relax Neck & Shoulders**-USB connector. The indicator light on the main unit lights red during charging, and purple when it is fully charged. The device is fully charged after 2 hours (2100 mAh).

Massage Function

EXPAIN Relax Neck & Shoulders consists of 4 innovative massage heads which can rotate in two directions. By rotating in one direction, the massage heads give soft massage. By rotating in the other direction the massage heads give hard massage.

1. Turn on the massage by holding the on / off button for 2 seconds.
2. Place the neck massager behind the neck or shoulders or where your muscle pain is.
3. The neck massager is now running a 15 minute massage program with active heating.
4. You can change the direction by a short press on the on/off button. The massage changes from hard to soft or soft to hard. Hard and soft massage is indicated by blue or purple light respectively.
5. If the massage stops, you may need to charge the battery or

- start the massage program again.
6. To stop the massage, hold the button for 2 seconds.
 7. **EXPAIN Relax Neck & Shoulders** has a flexible joint in the middle that allows the angle between the massage heads to be varied between 90° and 180°, therefore you can massage both the neck, shoulders, thighs and arms in many ways.

Heating function

Heat function is always on so **EXPAIN Relax Neck & Shoulders** provides for good blood circulation and is similar in feeling to a manual massage.

Battery and charging

Expain Relax Neck & Shoulders is recharged with the included USB charge cable and transformer.

The device has a battery life of about 1 hour during normal use. The device shows a low battery with red flashes. Use the unit for charging and charge it completely before removing the plug. You can use **EXPAIN Relax Neck & Shoulders** with a cord.

During charging **EXPAIN Relax Neck & Shoulders** lights red. When the device is fully charged the light is purple, and it is ready for cordless use again.

Maintenance and Storage

Clean the device with a dry or slightly damp cloth. The device is best kept in a sofa or comfortable chair. Make sure the device does not get wet. The device should be transported, stored and operated at a temperature between 0 ° C and 40 ° C, at a relative humidity between 30% and 85% and an atmospheric pressure between 500 and 1060 hPa. If the device is not working properly, stop all use and seek support on www.expain.com

Technical specifications	
Number of massage heads	4
Technology in use	Sensotouch Adaptive Massage, IR-heat
Power	2100mAh LiPo rechargeable battery or cable
Size (DxWxH)	11x24x10cm
Weight	0,8 kg
Safety standards	EN 60335, UN38.3, RoHS, REACH
Made in	China

Safety Information

Areas of use

For use on the neck, shoulders,
thighs and arms

Safety Precautions

Do not use this product if you have open wounds or very sensitive skin that cannot tolerate massage. Do not use this product if you have a back who cannot tolerate normal, manual massage.

Warning!

Do not use the device elsewhere on the body other than those places specified in this manual.

Expain AS
Karl Johans gate 13
0154 Oslo, Norway
www.expain.com

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