EXPAIN®

CHANGE NECK & SHOULDERS

USER MANUAL

- › English
- › Norsk
- › Dansk
- › Svensk
- › Suomi


EXPAIN®

Contents

- 1 Components
- 2 Introduction
- 3 Before use
- 4 How to mount the device
- 5 How to unmount the device
- 6 Usage
- 7 Usage with app
- 8 Biofeedback and habits
- 9 Battery and charging
- 10 Maintenance and storage
- 11 Technical specifications
- 12 Safety information

Components

- 1 Button
- 2 Light
- 3 USB connector
- 4 Fast snap connector
- 5 Adhesive patches snap
- 6 USB cable

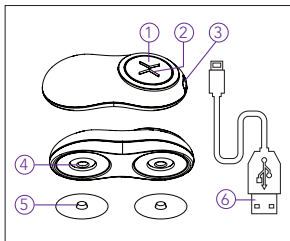Introduction

Congratulations on your new **EXPAIN Change Neck & Shoulders!** This is a unique device that allows you to monitor the muscle activity in the neck and shoulder muscles in a new way. **EXPAIN Change Neck & Shoulders** provides Biofeedback to help you change habits that may cause pain. Biofeedback is a learning process where observation of the signals from your own body, leads to a changed habits. The device vibrates when you have strained muscle activity patterns. This allows you to change your bad habits, which again may eliminate the cause of your pain.

EXPAIN Change Neck & Shoulders can be connected with a smartphone via Bluetooth Smart, or it can be used without connection to a phone. If you wish to keep track of your historical muscle activity pattern, download

EXPAIN app

(iOS or Android).

See separate section or www.expain.com for more information. The unit is wireless, powered with a rechargeable battery and is attached to the skin over the muscles in the neck / shoulder region with the included patch.

Before Use**EXPAIN Change Neck & Shoulders**

is made for personal use. Read instructions carefully. The unit must be charged before the first use. See section for battery and charging for more information on how the device is charged.


(NB: The device should not be used at night, and the unit should not be on the body continuously for over 16 hours. If there are signs of allergic reactions or rash, immediately stop using the device. Use the device with caution).

How to mount the device

1. Attach two patches to the bottom of the product with the metal snap buttons. Make sure both patches are securely connected to the product.
2. Remove the protective plastic film from both patches.
3. Mount the device on the body as indicated in the figure. Attach the device to the skin over the muscle on the side you have most pain. For most people this is the dominant side (Right for right-handers and left for left-handers). Place the product with the USB connector pointing downwards towards the arm. If you

wish to measure
on both sides
simultaneously,
we recommend that you
use one extra device.

Let the device be on
your body for the
time period you
want biofeedback.

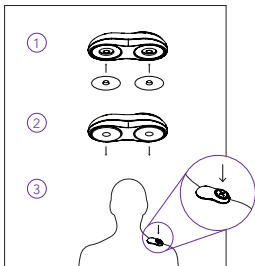How to unmount your device


1. Loosen the unit from the body by gently pulling the patches up from the skin. If the patches have lost much of their adhesive effect, they can be washed in water without soap. **NB!** *Remember to take patches off the **EXPAIN Change Neck & Shoulders** device before you wash. The unit is not waterproof.*
2. Recharge the unit to ensure that it is fully charged for the next use. Place the device back in the travel case after use for safe storage.

Usage

EXPAIN Change Neck & Shoulders

is turned on by pressing the button on top of the unit.

When using the device, it will vibrate when it recognizes strained muscle activity patterns. Read more in the section of Biofeedback and habits.

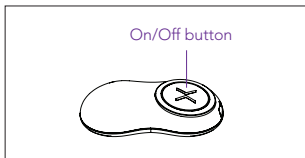Usage with app

You can get more from the **EXPAIN Change Neck & Shoulders** device if you connect it to a smartphone.

The unit will then automatically transfer data to your smartphone when it is within reach and you can follow the historical development in your app. You can, for example do ergonomic changes by looking at your muscle activity.

1. **EXPAIN Change Neck & Shoulders** are connectable whenever the device is turned on. It is not possible to connect the device to multiple phones simultaneously.

If you want to pair the device with another phone, you must remove the device in the settings in your

app / phone, and restore factory settings if necessary, see the Battery and charging section.

The **EXPAIN app** can be downloaded from Apple App Store or Google Play Store. Follow the instructions in the app. The app has a built in pain diary that allows you to compare your muscle activity patterns with your pain. In this way you can see when you've had the most strained muscle activity and if this has a relationship with your pain. It is recommended to register your pain regularly, and to use the device at least for a few days to get the most out of **EXPAIN Change Neck & Shoulders**.

Light and vibration

The table below shows the meaning of various lights and vibrations.

Function	Light/Vibration
Connected Bluetooth Smart	Blue flash
Low battery	Red flash
On, without Bluetooth connection	Purple flash
Low battery	Red flash
Fully charged	Purple
Notice sustained muscle activity	Short vibrations
Low battery	Long vibrations

Biofeedback and habits

Habit change happens when the device is used often and you take breaks or do some exercises every time the device vibrates. Research shows that biofeedback, which is a learning process, helps to improve health and overall performance by observing the signals generated by one's own body.

The device will vibrate when it recognizes strained muscle activity patterns (Sustained muscle activity). When the device vibrates, try to take a short break from the work, do some shoulder exercises or lower your shoulders so that the muscles in neck / shoulder regions relax. A good relaxation exercise is to breathe calmly in with the nose and then releasing the breath out of

your mouth while you let your arms and shoulders hang down by the side of the body. Repeat this exercise 2-5 times.

Battery and charging

The device has a battery life of approximately 8 hours of normal use. The sensor flashes red to indicate low battery. In addition, the device vibrates Long


vibrations when the battery is critically low. **EXPAIN Change**

Neck & Shoulders is charged with supplied or corresponding micro

USB cable. Unmount the device from your

body and charge it completely before next use.

The light is red during charging, and purple when it is fully charged.

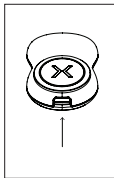Restoring factory settings

If you can't connect your phone the device can be reset to its original state. This is done as follows: Turn the device on and set it to charge. Press the button until the device starts flashing red. When the device lights red and vibrates, the factory settings are restored.

Maintenance and Storage

Clean the device with a dry or slightly damp cloth. The device is best kept in the accompanying travel case. Also carry your extra patches and USB charging cable in this case so that you always have them available.

Make sure the device does not get wet. It is not recommended to use the device in humid environments or if you plan to do activities that will

cause sweating.

The patches can be worn out, but are made to last at least 10 times / 10 hours of use. If the adhesive effect is low, the patches can be gently washed in water. **NB!** *Remember to remove the patches from the device before washing, otherwise the unit can be destroyed by moisture penetration. NB! If the desired adhesive effect is not obtained by washing, the patches are ready for replacement.*

EXPAIN Change Neck & Shoulders comes with three pairs of patches. New patches can be purchased at www.expain.com

The unit should be transported, stored and operated in temperatures between 10°C and 40°C, at a relative humidity

between 30% and 85% and with an atmospheric pressure between 500 and 1060 hPa.

Keep the device and patches in their original packaging or supplied travel case when not in use.

Cover patches with the protective film when they are not in use. If the device does not work correctly, contact support at www.expain.com

Tehnickal specifications	
Technology	Surface Electromyography (sEMG)
Power	LiPo rechargeable battery
Bluetooth	Bluetooth 4.2 Smart
Size (DxBxH)	32 mm x 60 mm x 13 mm
Weight	20 g
Safety	IEC 62133 EN 60950
EMC	EN 301 489-17 EN 300 328 EN 301 489-1
Safety	RoHS, REACH Primary Skin Irritation Test
Made in	Norway

Safety information

Applications

For Biofeedback against strained (sustained) muscle activity in the Upper Trapezius muscle.

Safety precautions

Do not use this product if you have very sensitive skin.

Warning!

Do not use the device on other locations on the body other than those specified in this manual.

Expain AS
Karl Johans gate 13
0154 Oslo, Norway
www.expain.com

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